


Amino-Modifier-15-dT CEP
Product No. BA 0289
Product Information

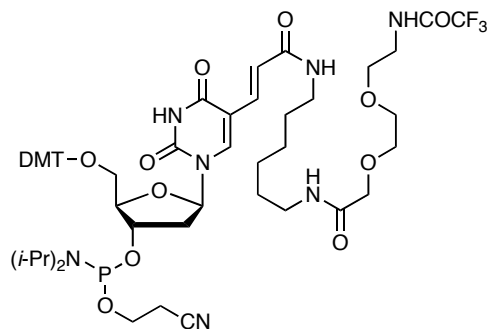

$C_{56}H_{73}F_3N_7O_{13}P$
Mol. Wt.: 1140.19

Amino-Modifier-15-dT CEP (BA 0289) is similar to Amino-Modifier-C6-dT CEP (BA 0015) except that it offers a longer tether that includes an amphipathic glycol ether region.

Use: Employ acetonitrile diluent at the concentration recommended by the synthesizer manufacturer. Use standard coupling protocols; extended coupling is not required. Cleavage from the solid support and nucleobase deprotection with concentrated ammonium hydroxide may be carried out using standard protocols.